

Re-Order #	Description	Size	U/M
703286	Econo Footrest	16"-18" x 1"	Ea.
703288	Econo Footrest	16"-18" x 2"	Ea.
703294	Econo Footrest - Bariatric	20"-24" x 1"	Ea.
703296	Econo Footrest - Bariatric	20"-24" x 2"	Ea.
703297	Econo Footrest - Bariatric	20"-24" x 3"	Ea.

Note: All models above fit 16"-18" and 20"-24" Wheelchairs and Geri-Chairs



Purpose

Skil-Care's Econo Footrest Extender is designed to solve two important problems: It provides support to control foot drop and it prevents residents' feet from slipping off the footrests. The Econo Footrest Extender is designed to fit 16-inch, 18-inch and 20-inch wide wheelchairs and geri-chairs depending on the model chosen.

Application

1. Before seating resident, lock wheel brakes and rotate footrests to the out-of-the-way position. Assist patient into chair and position him/her as far back on the seat as possible. Reposition footrests.
2. Place Econo Footrest Extender on footrests. Cushion section on footrest, padded backboard with straps against footrest support bars (*Fig. 1*).
3. Place straps behind footrest support bars, secure buckles and tighten straps (*Fig. 2*).
4. The Econo Footrest Extender may be positioned to accommodate various degrees of ankle plantar flexion. This may be accomplished by securing the top straps to the upper section of the footrest support bars as shown in the illustration (*Fig. 3*).
5. Important Safety Note. Remove the Econo Footrest Extender and rotate the footrests to the out-of-the-way position when assisting your resident from the chair. Do not permit your resident to stand with footrests and Econo Footrest Extender in place.

Easy Care

Skil-Care's Econo Footrest Extender is made from durable vinyl for longlasting use. To clean, just spray and wipe with a damp cloth.

Fig. 1

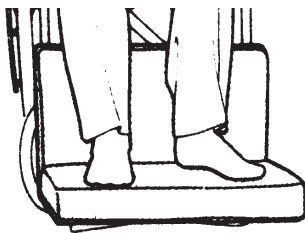


Fig. 2

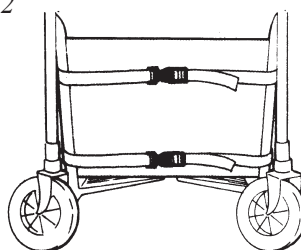


Fig. 3

