
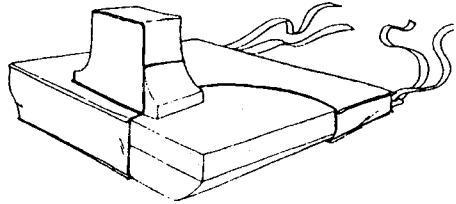


# LEG ABDUCTOR

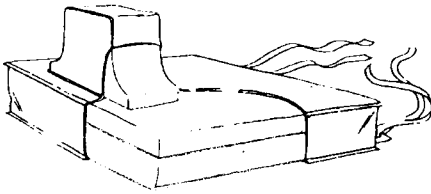

**LEVEL - BASE LEG ABDUCTOR**

**WEDGE - BASE LEG ABDUCTOR**

Reorder #	Description	Cover	Size (Inches)
☐758010	Level-Base Abductor/Foam	Vinyl	16" x 18" x 3½"
☐758015	Level-Base Leg Abductor/Gel-Foam	Vinyl w/Cloth	16" x 18" x 3½"
☐758020	Wedge-Base Leg Abductor/Foam	Vinyl	16" x 18" x 4½ x 1½"
☐758025	Wedge-Base Leg Abductor/Gel-Foam	Vinyl w/Cloth	16" x 18" x 4½ x 1½"

The contoured foam pommel keeps legs in abduction without compromising comfort. The firm foundation eliminates sling-seat hammocking to control inner thigh rotation, improve sitting posture, and reduce pressure on the spine. Available in a foam-only version for basic comfort and pressure reduction, and a gel-foam version when additional pressure relief is required. Use the level-base model for most patients, or the wedge-base model when forward sliding is a problem.

## Important Safety Feature

All Skil-Care Leg Abductors have safety ties for the purpose of securing the abductors to the wheelchair. These ties must be used to prevent the abductors from slipping off wheelchair seats.

## Cover Care

The incontinent-proof vinyl cover can be cleaned with a spray-and-wipe cleanser and damp cloth. A mild disinfectant can be used. Gel-foam models have outer cloth covers that can be laundered and dried at temperatures that do not exceed 180°F (82°C).

Do not use bleach.