

The **Premier Pressure-Check dual-layer, 3-zone 6" Mattress** is specifically designed to provide an industry leading pressure redistribution and pressure relief environment for the resident. The top layer is made out of a unique polyurethane foam material and is computer-cut into individual "Support Pods". The "Support Pod" technology optimally addresses the key issues of weight distribution, pressure relief, anti-shearing, and provides a positive thermoclimate by eliminating heat buildup. The 400 individual "Support Pods" separate the mattress into three specific pressure redistribution zones; Head, Heel and Trunk. The bottom layer is high-density compression-resistant foam to support weight distribution and extend the mattress life. The Low-Shear III cover is a 4-Way Stretch low moisture, super-durable, breathable, fluid-proof fabric that wipes clean for easy care while protecting sensitive skin against friction and shear damage. The bottom section of the mattress cover is slip-resistant for the patients safety.

## Pressure-Check Mattress

**Pressure-Check™ Foam** has 400 support pods that adjust individually for optimal weight distribution and effective pressure relief.


**Heel Zone** features smaller foam sections designed to accommodate the special pressure relief needs of sensitive heels.

**High Density Foam Base** provides uniform, compression-resistant support for weight distribution, comfort and extended mattress life.

**Waterfall Zipper** prevents liquid penetration


The **NEW Low-Shear III** cover provides an industry top choice of 4-Way Stretch fabric that provides less moisture, more comfort, and cooler feel than a standard 2-Way stretch fabric.

LSIII Four-Way Stretch

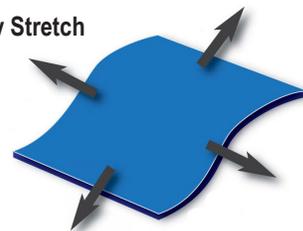


## Gel-Infused Visco Mattress

For residents who want additional comfort we recommend the gel-infused visco mattress. The biggest advantage of the Gel-infused mattress is that it fully cradles to the body and keeps the weight off of the pressure points while the gel in the foam provides a cooler surface by allowing the body heat to escape.

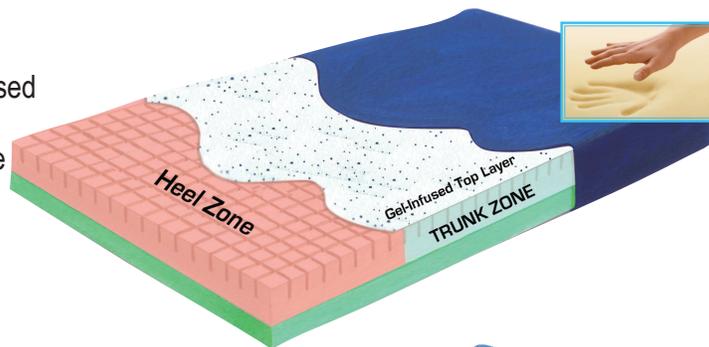

## Perimeter Guard Mattress

Mattresses are sold with or without perimeter guard. The perimeter guard mattresses come with a two inch (2") raised perimeter and provides a safer transfer support and act as a gentle reminder that the resident is near the mattress edge. Skil-Care's EZ exit cutout prevents user injury and is safer for the caregiver.

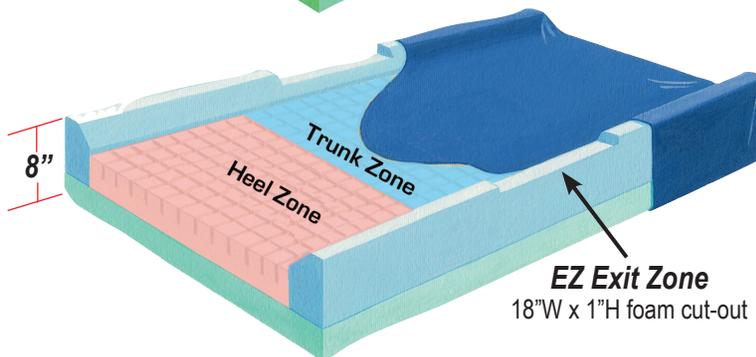

**Whatever mattress model you choose rest-assured, Skil-Care has you covered!**

All Skil-Care multi-layer mattresses deliver a full six inches (6") of pressure redistribution with enhanced scapular, sacral and heel section to further reduce pressure in three vulnerable areas as shown in the pressure mapping chart.


All of the mattresses comply with fire code 16 CFR 1633 and are quality built to last.

All mattresses come with a Seven (7) year non-prorated replacement warranty.

Foam	Gel-Infused Visco	Dimensions	Weight (max.)
558015	558123	6" x 36" x 76"	500 lbs.
558020	558125	6" x 36" x 80"	500 lbs.
558025	558127	6" x 36" x 84"	500 lbs.

**Perimeter Guard**

Foam	Gel-Infused Visco	Dimensions	Weight (max.)
558051	558133	6" x 36" x 76"	500 lbs.
558052	558135	6" x 36" x 80"	500 lbs.
558053	558137	6" x 36" x 84"	500 lbs.

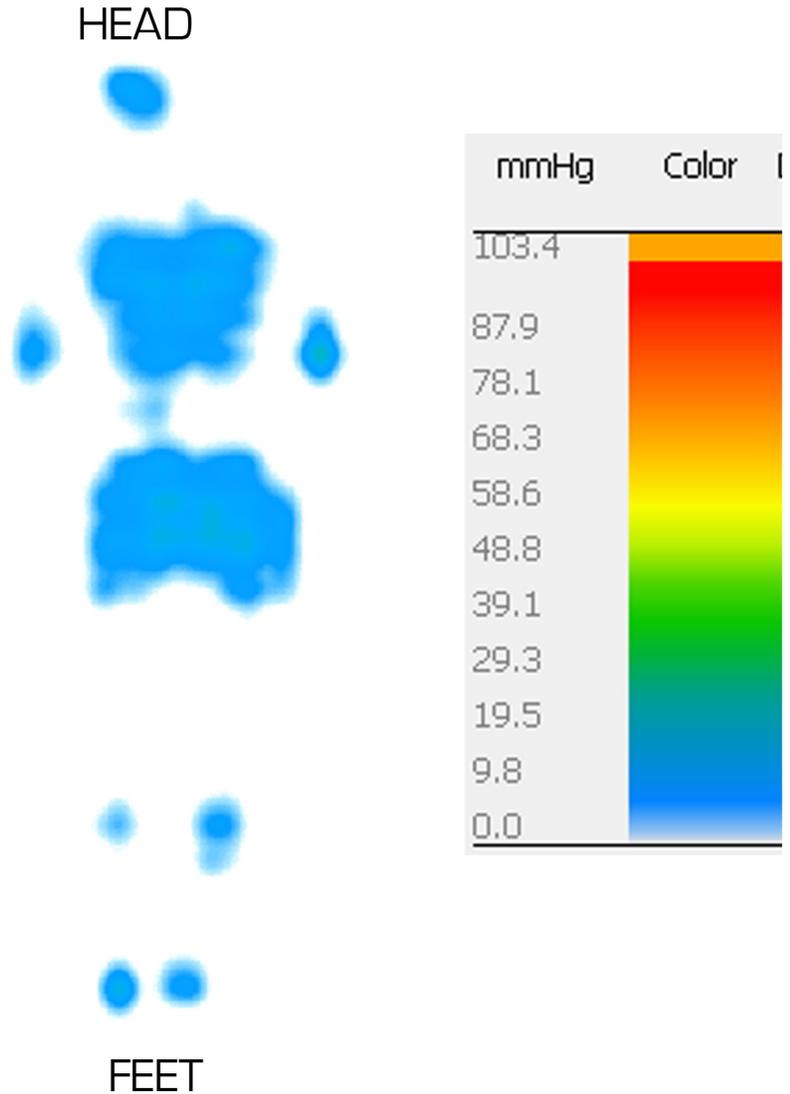

Skil-Care custom manufacturers Bariatric mattresses to your specifications and offers Psychiatric mattresses with a **NO ZIPPER** cover.

**For any questions, information or pricing, please contact our customer service department at 800.431.2972**

**Visit us online [www.skil-care.com](http://www.skil-care.com)**