
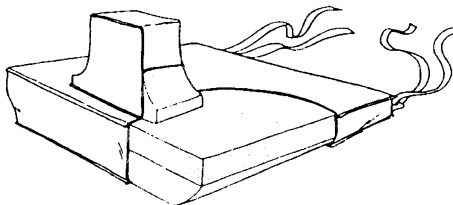


LEG ABDUCTOR

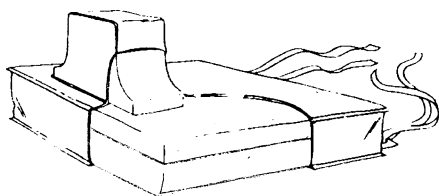

LEVEL - BASE LEG ABDUCTOR

WEDGE - BASE LEG ABDUCTOR

Reorder #	Description	Cover	Size (Inches)
☐758010	Level-Base Abductor/Foam	Vinyl	18" x 16" x 3½"
☐758015	Level-Base Leg Abductor/Gel-Foam	Vinyl w/LSI	18" x 16" x 4"
☐758020	Wedge-Base Leg Abductor/Foam	Vinyl	18" x 16" x 4" x 1½"
☐758025	Wedge-Base Leg Abductor/Gel-Foam	Vinyl w/LSI	18" x 16" x 4½ x 2"

The contoured foam pommel keeps legs in abduction without compromising comfort. The firm foundation eliminates sling-seat hammocking to control inner thigh rotation, improve sitting posture, and reduce pressure on the spine. Available in a foam-only version for basic comfort and pressure reduction, and a gel-foam version when additional pressure relief is required. Use the level-base model for most patients, or the wedge-base model when forward sliding is a problem.

Important Safety Feature

Skil-Care Leg Abductors have safety ties for the purpose of securing the abductors to the wheelchair attached to the vinyl inner cover. If you choose the LSI cover model number, you will have a friendly release buckle attachment. These attachments must be used to prevent the abductors from slipping off wheelchair seats.

Cover Care

The incontinent-proof vinyl cover can be cleaned with a spray-and-wipe cleanser and damp cloth. A mild disinfectant can be used. Gel-foam models have outer cloth covers that can be laundered and dried at temperatures that do not exceed 180°F (82°C). Do not use bleach.